



The professional cell culture empowers a healthier world

# **DMEM (High glucose) Complete Medium (10)**

Cat. No.: PM150210B

Size: 125mL×4

### **General Information**

**Product Form** Liquid

Concentration Ready-to-use

Components DMEM (High glucose)[PM150210]+10% Nutrient+1% Supplement

**Bacterial detection** Negative Fungal detection Negative Mycoplasmal detection Negative **Endotoxin level** < 3 EU/mL Shipping Ice bag

2-8°C, Shading Light Storage

**Expiry date** 6 months

## **Product Introduction**

icella DMEM (Dulbecco's Modified Eagle Medium) was developed on the basis of MEM medium. Compared with MEM medium, the content of amino acid increased by 2 times, the content of vitamin increased by 4 times, and the content of non-essential amino acid, trace iron ion and sodium pyruvate were increased by 4 times.

The glucose content of DMEM medium was originally designed as 1000 mg/L (low Glucose type), and then developed into 4500 mg/L (high Glucose type), which has been widely used in cell culture.

DMEM (High glucose) was widely used in fast growth, low adhesion cells, hybridoma myeloma cells, clone cells, DNA transfected transformation cells, various primary virus host cells, single cell culture and vaccine production.

#### **Guidelines for use**

- 1. Pricella's cell culture media undergoes strictguality control to ensure sterility, but may get contaminated during use. Follow these guidelines for sterile handling toavoid contamination.
- 2. Always wipe your gloved hands and work area with 70% ethanol.
- 3. Wipe the outside of the containers, flasks, plates, and dishes with 70% ethanol before placing them in the cell culture hood.
- 4. Use sterile pipette tips and pipettes to work with liquids, and use each pipette tip only once to avoid cross contamination. Do not unwrap sterile pipettes until they are ready to beused. Keep pipettes and tips within the clean work area.
- 5. Do not talk while performing sterile procedures and perform your cell culture as efficiently and carefully as possible tominimize contamination.

# **Quality control**

Standard evaluations for cell culture media are pH, osmolality, endotoxins and sterility testing for liquid products, cell growth experiments.

### **Notes**


Tel:1-832-243-6086 Fax:1-832-243-6017 Toll-free:1-888-852-8623 Email: techsupport@pri-cella.com Web: www.pri-cella.com

Rev. V3.0





- 1. This product is for research use only.
- 2. This product is sterilized by 0.1  $\mu m$  filtration.
- 3. It is necessary to pay attention to the aseptic operation and avoid the contamination during the culture.
- 4. It is not suitable for long time storage at room temperature.
- 5. This product is a ready-to-use medium. If there is no special need, don't add serum, penicillin and streptomycin. It can be used directly.

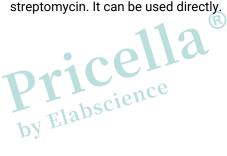







Tel:1-832-243-6086 Fax:1-832-243-6017 Toll-free:1-888-852-8623 Web: www.pri-cella.com