Elabscience®


Elab Fluor[®] 647 Anti-Mouse CD49b/pan-NK cells Antibody[DX5]

Catalog Number: E-AB-F1116M

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description	
Reactivity	Mouse
Host	Rat
lsotype	Rat IgM, κ
Clone No.	DX5
Isotype Control	Elab Fluor [®] 647 Rat IgM, κ Isotype Control[RTK2118] [Product E-AB-F09772M]
Conjugation	Elab Fluor [®] 647
Conjugation Information	Elab Fluor [®] 647 is designed to be excited by the Red laser (627-640 nm) and detected using an optical filter centered near 670 nm (e.g., a 660/20 nm bandpass filter).
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Applications	Recommended usage
FCM	Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 µL of antibody per test (million cells in 100 µL staining volume or per 100 µL of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Data

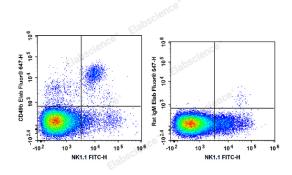

C57BL/6 murine splenocytes are stained with FITC Anti-

Mouse CD161/NK1.1 Antibody and Elab Fluor[®] 647 Anti-Mouse CD49b Antibody (Left). Splenocytes are stained with

FITC Anti-Mouse CD161/NK1.1 Antibody and Elab Fluor[®] 647 Rat IgM, κ Isotype Control (Right).

Preparation & Storag	je
Storage	Keep as concentrated solution.
	This product can be stored at 2-8°C for 12 months. Please protected from prolonged exposure to light and do not freeze.
Shipping	Ice bag
Antigen Information	
Alternate Names	CD49 antigen-like family member B;CD49b;Collagen receptor;GPla;Integrin alpha-2;
	Platelet membrane glycoprotein la;VLA-2 subunit alpha;pan-NK cells

For Research Use Only

Toll-free: 1-888-852-8623 Web:<u>w w w .elabscience.com</u> Tel: 1-832-243-6086 Email:techsupport@elabscience.com

Elabscience®

Uniprot ID Gene ID Background Q62469 16398

DX5 antigen has been recently characterized as CD49b. It is a 150 kD integrin α chain also known as α 2 integrin, VLA-2 α chain, and integrin α 2 chain. CD49b non-covalently associates with CD29 (β 1 integrin) to form the CD49b/CD29 complex known as VLA-2, a receptor for collagen and laminin. CD49b is expressed on platelets, the majority of NK cells, NKT cells, and a small subset of CD8+ T cells (this population can be significantly increased following viral infection). DX5 is used for the identification and isolation of NK cells, and is especially useful for identifying NK cells in mice lacking the NK1.1 antigen.