Recombinant Mouse P-selectin/CD62P Protein (Fc Tag)

Catalog Number: PKSM041284



Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description	
Species	Mouse
Mol_Mass	99.5 kDa
Accession	Q01102

Bio-activity Not validated for activity

т.				
P	rn	ne	rti	AC
ж.	LU			

Purity > 95 % as determined by reducing SDS-PAGE.

Endotoxin < 1.0 EU per µg of the protein as determined by the LAL method.

Storage Generally, lyophilized proteins are stable for up to 12 months when stored at -20 to -80

°C. Reconstituted protein solution can be stored at 4-8°C for 2-7 days. Aliquots of

reconstituted samples are stable at < -20°C for 3 months.

Shipping This product is provided as lyophilized powder which is shipped with ice packs.

Formulation Lyophilized from a 0.2 μm filtered solution of PBS, pH 7.4.


Normally 5% - 8% trehalose, mannitol and 0.01% Tween 80 are added as protectants

before lyophilization.

Please refer to the specific buffer information in the printed manual.

Reconstitution Please refer to the printed manual for detailed information.

Data

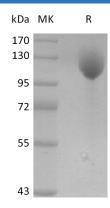

Background

P-selectin/CD62P is a single-pass type I membrane protein which is a member of the Selectin family. It consists 768 amino acid (aa). P-selectin is a cell surface glycoprotein expressed by activated platelets and endothelial cells. It induced expression in lung, liver, kidney and heart after endotoxin treatment. Ca2+-dependent receptor for myeloid cells that binds to carbohydrates on neutrophils and monocytes. It mediates the interaction of activated endothelial cells or platelets with leukocytes. The ligand recognized is sialyl-Lewis X. it also mediates rapid rolling of leukocyte rolling over vascular surfaces during the initial steps in inflammation through interaction with PSGL1. P-selectin interacts with SNX1 7, PSGL1/SEPL, PODXL2, mediates neutrophil adhesion and leukocyte rolling. This interaction requires the sialyl-Lewis X epitope of PSGL1 and PODXL2, and specific tyrosine sulfation on PSGL1.

For Research Use Only