
FITC Anti-Human/Mouse/Rat CD47 Antibody[MIAP410]

Catalog Number: E-AB-F1016C

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description	
Reactivity	Human;Mouse;Rat
Host	Mouse
Isotype	Mouse IgG1, ĸ
Clone No.	MIAP410
Isotype Control	FITC Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792C]
Conjugation	FITC
Conjugation Information	FITC is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 530 nm (e.g., a 525/40 nm bandpass filter).
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Applications	Recommended usage
FCM	Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Data

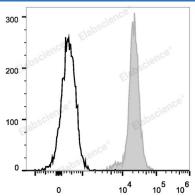

C57BL/6 murine splenocytes are stained with FITC Anti-Human/Mouse/Rat CD47 Antibody (filled gray histogram). Unstained splenocytes (empty black histogram) are used as control.

Preparation & Storag	le
Storage	Keep as concentrated solution.
	This product can be stored at 2-8°C for 12 months. Please protected from prolonged
	exposure to light and do not freeze.
Shipping	Ice bag
Antigen Information	
Alternate Names	Cd47;IAP;Integrin-associated protein;Leukocyte surface antigen CD47
Uniprot ID	Q08722;Q61735;P97829
Gene ID	961.16423.29364

For Research Use Only

Toll-free: 1-888-852-8623 Web:<u>w w w .elabscience.com</u>

Elabscience®

Elabscience Bionovation Inc. A Reliable Research Partner in Life Science and Medicine

Background

CD47, also known as Integrin-Associated Protein (IAP), is a membrane protein of about 50 kD with an IgV-like extracelluluar domain, a five membrane-spanning segment and a short terminal cytoplasmic region. It is widely expressed on many cell types and often associated with beta 3 integrins. It has been reported that CD47 functions as a self marker. Red cells lacking CD47 were rapidly cleared from the bloodstream by splenic macrophages. By binding to SIRPα, CD47 controls hemostatic innate immune functions, such as phagocytosis and cell trafficking.

For Research Use Only