

Elab Fluor® 700 Anti-Mouse CD49b/pan-NK cells Antibody[DX5]

Catalog Number: E-AB-F1116UM1

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description

Reactivity	Mouse
Host	Rat
Isotype	Rat IgM, κ
Clone No.	DX5
Isotype Control	Elab Fluor® 700 Rat IgM, κ Isotype Control[RTK2118] [Product E-AB-F09773M1]
Conjugation	Elab Fluor® 700
Conjugation Information	Elab Fluor® 700 is designed to be excited by the Red laser (627-640 nm) and detected using an optical filter centered near 719 nm (e.g., a 725/40 nm bandpass filter).
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% sodium azide and 1% BSA.


Applications

Recommended usage

FCM

Each lot of this antibody is quality control tested by flow cytometric analysis. Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use. We suggest each investigator should titrate the reagent to obtain optimal results [The recommended concentration is 0.1-1 μg/10⁶ cells in 100 μL volume].

Data

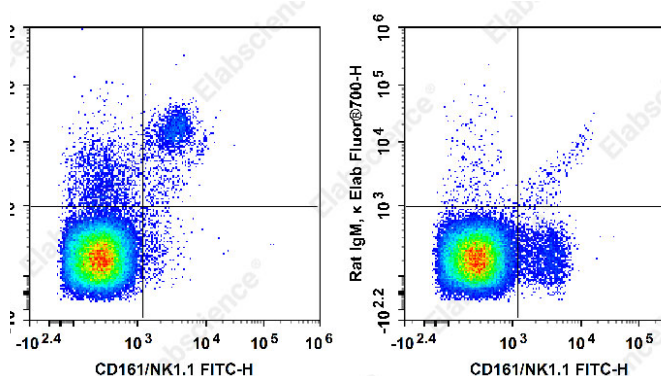

Staining of C57BL/6 murine splenocytes with FITC Anti-

Mouse CD161/NK1.1 Antibody[PK136] and Elab Fluor® 700 Anti-Mouse CD49b/pan-NK cells Antibody[DX5](left) or Elab

Fluor® 700 Rat IgM, κ Isotype Control(right). Total viable cells were used for analysis.

Preparation & Storage

Storage	Keep as concentrated solution. This product can be stored at 2-8°C for 12 months. Please protected from prolonged exposure to light and do not freeze.
Shipping	Ice bag

Antigen Information

Alternate Names	CD49 antigen-like family member B;CD49b;Collagen receptor;GPIa;Integrin alpha-2; Platelet membrane glycoprotein Ia;VLA-2 subunit alpha;pan-NK cells
Uniprot ID	Q62469

For Research Use Only

Gene ID

16398

Background

DX5 antigen has been recently characterized as CD49b. It is a 150 kD integrin α chain also known as $\alpha 2$ integrin, VLA-2 α chain, and integrin $\alpha 2$ chain. CD49b non-covalently associates with CD29 ($\beta 1$ integrin) to form the CD49b/CD29 complex known as VLA-2, a receptor for collagen and laminin. CD49b is expressed on platelets, the majority of NK cells, NKT cells, and a small subset of CD8+ T cells (this population can be significantly increased following viral infection). DX5 is used for the identification and isolation of NK cells, and is especially useful for identifying NK cells in mice lacking the NK1.1 antigen.