

A Reliable Research Partner in Life Science and Medicine

PE Anti-Mouse TER-119 Antibody[TER-119]

Catalog Number: E-AB-F1125UD

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description

Reactivity Mouse Host Rat

Isotype Rat IgG2b, κ **Clone No.** TER-119

Isotype Control PE Rat IgG2b, κ Isotype Control[LTF-2] [Product E-AB-F09843D]

Conjugation PE

Conjugation Information PE is designed to be excited by the Blue (488 nm), Green (532 nm) and Yellow-Green

(561 nm) lasers and detected using an optical filter centered near 575 nm (e.g., a 585/42

nm bandpass filter).

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer.

Applications Recommended usage


FCM Each lot of this antibody is quality control tested by flow cytometric analysis. Please

check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use. We suggest each investigator should titrate the

reagent to obtain optimal results [The recommended concentration is 0.1-1 μg/10⁶ cells

in 100 µL volume].

Data

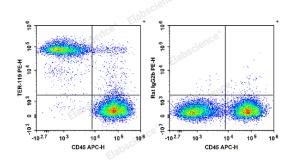

C57BL/6 murine bone marrow cells are stained with APC Anti-Mouse CD45 Antibody and PE Anti-Mouse TER-119 Antibody (Left). Bone marrow cells are stained with APC Anti-Mouse CD45 Antibody and PE Rat IgG2b, κ Isotype Control (Right).

Preparation & Storage

Storage Keep as concentrated solution.

This product can be stored at 2-8°C for 12 months. Please protected from prolonged

exposure to light and do not freeze.

Shipping Ice bag

Antigen Information

Alternate Names Ly-76;Lymphocyte antigen 76;TER119

Gene ID 104231

For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086 Fax: 1-832-243-6017

Web:www.elabscience.com Email:techsupport@elabscience.com

Elabscience Bionovation Inc.

A Reliable Research Partner in Life Science and Medicine

Background

Elabscience®

The TER-119 antigen is a 52 kD glycophorin A-associated protein, also known as Ly-76. TER-119 is an erythroid-specific antigen expressed on early proerythroblasts to mature erythrocytes, but not on erythroid colony-forming cells (BFU-E, blast-forming unit erythroid, or CFU-E, colony-forming unit erythroid).

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086
 Fax: 1-832-243-6017

 Web:www.elabscience.com
 Email:techsupport@elabscience.com