

Elab Fluor® 647 Anti-Mouse Ig light chain κ Antibody[RMK-45]

Catalog Number: AN00649M

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description


Reactivity	Mouse
Host	Rat
Isotype	Rat IgG
Clone No.	RMK-45
Isotype Control	Elab Fluor® 647 Rat IgG1, κ Isotype Control[HRPN] [Product E-AB-F09822M]
Conjugation	Elab Fluor® 647
Conjugation Information	Elab Fluor® 647 is designed to be excited by the Red laser (627-640 nm) and detected using an optical filter centered near 670 nm (e.g., a 660/20 nm bandpass filter).
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer.

Applications

Recommended usage

FCM	Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μL of antibody per test (million cells in 100 μL staining volume or per 100 μL of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.
------------	---

Data

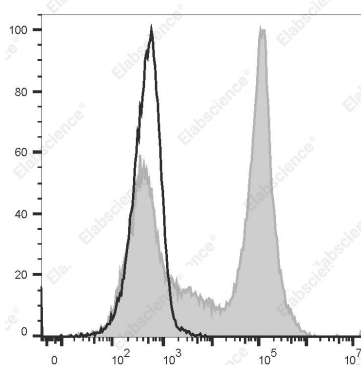

Staining of C57BL/6 murine splenocytes cells with Elab Fluor

® 647 Anti-Mouse Ig light chain κ Antibody[RMK-45] (filled gray histogram) or Elab Fluor® 647 Rat IgG Isotype Control (empty black histogram). Total viable cells were used for analysis.

Preparation & Storage

Storage	Keep as concentrated solution. This product can be stored at 2-8°C for 24 months. Please protected from prolonged exposure to light and do not freeze.
Shipping	Ice bag

Antigen Information

Alternate Names	Immunoglobulin light chain kappa
Gene ID	243469

For Research Use Only

Toll-free: 1-888-852-8623
Web: www.elabscience.com

Tel: 1-832-243-6086
Email: techsupport@elabscience.com

Fax: 1-832-243-6017

Rev. V1.2

Background

The RMK-45 monoclonal antibody reacts with immunoglobulin light chain kappa in all tested mouse haplotype (Igh-a and b). It does not react with the lambda chain.

For Research Use Only

Toll-free: 1-888-852-8623
Web: www.elabscience.com

Tel: 1-832-243-6086
Email: techsupport@elabscience.com

Fax: 1-832-243-6017

Rev. V1.2