A Reliable Research Partner in Life Science and Medicine

Recombinant Rat E-Selectin/SELE Protein (Fc Tag)

Catalog Number: PKSR030388

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description

Species Rat

Source HEK293 Cells-derived Rat E-Selectin/SELE protein Met1-Pro494, with an C-terminal

Calculated MW 78.8 kDa Observed MW 114-119 kDa Accession NP 620234.1

Measured by the ability of the immobilized protein to support the adhesion of U937 **Bio-activity**

> cells. When 5 x 10^4 cells/well are added to SELE-coated plates (2 μ g/mL and 100 μ L/well), approximately > 85% cells will adhere specifically after 60 minutes at 37°C.

Properties

Purity > 95 % as determined by reducing SDS-PAGE.

Endotoxin < 1.0 EU per µg of the protein as determined by the LAL method.

Generally, lyophilized proteins are stable for up to 12 months when stored at -20 to -80 Storage

°C. Reconstituted protein solution can be stored at 4-8°C for 2-7 days. Aliquots of

reconstituted samples are stable at < -20°C for 3 months.

This product is provided as lyophilized powder which is shipped with ice packs. Shipping

Lyophilized from sterile PBS, pH 7.4 Formulation


Normally 5% - 8% trehalose, mannitol and 0.01% Tween 80 are added as protectants

before lyophilization.

Please refer to the specific buffer information in the printed manual.

Reconstitution Please refer to the printed manual for detailed information.

Data

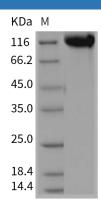

> 95 % as determined by reducing SDS-PAGE.

Background

Elabscience Bionovation Inc.



A Reliable Research Partner in Life Science and Medicine

E-selectin, also known as endothelial leukocyte adhesion molecule-1 (ELAM-1) and CD62E, is an inducible adhesion molecule that is expressed on the surfaces of stimulated vascular endothelial cells and is sometimes involved in cancer cell metastasis. As a member of the LEC-CAM or selectin family, E-selectin recognises and binds to sialylated carbohydrates including members of the Lewis X and Lewis A families found on monocytes, granulocytes, and T-lymphocytes. E-selectin supports rolling and stable arrest of leukocytes on activated vascular endothelium, and furthermore, it was indicated that it can also transduce an activating stimulus via the MAPK cascade into the endothelial cell during leukocyte adhesion. E-selectin regulates adhesive interactions between certain blood cells and endothelium. The soluble form of E selectin (sE-selectin) is a marker of endothelial activation, and has a potential role in the pathogenesis of cardiovascular disease as raised levels have been found in hypertension, diabetes and hyperlipidemia, although its association in established atherosclerosis disease and its value as a prognostic factor is more controversial. soluble E-selectin is inversely associated with the muscular component of the left ventricle, thereby suggesting that the lack of such a reparative factor may be associated with cardiac remodeling in end-stage renal disease (ESRD) patients. In addition, this adhesion molecule appears to be involved in the pathogenesis of atherosclerosis.

Fax: 1-832-243-6017