Elabscience®

PE Anti-Mouse CD51 Antibody[RMV-7]


Catalog Number: E-AB-F1235UD

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description	
Reactivity	Mouse
Host	Rat
lsotype	Rat lgG1, κ
Clone No.	RMV-7
Isotype Control	PE Rat IgG1, κ Isotype Control[HRPN] [Product E-AB-F09823D]
Conjugation	PE
Conjugation Information	PE is designed to be excited by the Blue (488 nm), Green (532 nm) and Yellow-Green (561 nm) lasers and detected using an optical filter centered near 575 nm (e.g., a 585/42 nm bandpass filter).
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% sodium azide and 1% BSA.
Applications	Recommended usage
FCM	Each lot of this antibody is quality control tested by flow cytometric analysis. Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use. We suggest each investigator should titrate the

reagent to obtain optimal results [The recommended concentration is 0.1-1 μ g/10⁶ cells in 100 μ L volume].

Data

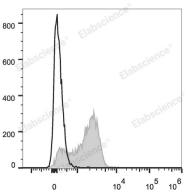

C57BL/6 murine bone marrow cells are stained with PE Anti-Mouse CD51 Antibody (filled gray histogram). Unstained bone marrow cells (empty black histogram) are used as control.

Preparation & Storage	
Storage	Keep as concentrated solution.
	This product can be stored at 2-8°C for 12 months. Please protected from prolonged exposure to light and do not freeze.
Shipping	Ice bag
Antigen Information	
Alternate Names	ITGAV;Integrin alpha-V;Integrin αV chain;Vitronectin Receptor;αV integrin
Uniprot ID	P43406
Gene ID	16410

For Research Use Only

Elabscience®

Background

CD51 is a 140 kD protein, also known as α V integrin, vitronectin receptor, and integrin α V. It is a member of the integrin family, expressed on activated T cells, polymorphonuclear granulocytes, platelets, blastocysts, and osteoclasts. CD51 forms heterodimers by association with integrins β 1, β 3, β 5 or β 6; these complexes then act as receptors for multiple extracellular matrix proteins (ECM). The α vintegrin heterodimers have varied functions in development, stimulation/activation and homeostasis. The primary ligands for CD51 complexes are fibronectin, fibrinogen, vitronectin, thrombspondin, von Willebrand factor, and CD31. The RMV-7 antibody has been reported to block binding of CD51 to vitronectin, fibronectin, and CD31 in some cell types, as well as blocking LAK cell cytotoxicity.