A Reliable Research Partner in Life Science and Medicine

# Recombinant Human FGF-2/FGF basic/FGFb Protein

Catalog Number: PKSH033503

Note: Centrifuge before opening to ensure complete recovery of vial contents.

### **Description**

Species Human

**Source** E.coli-derived Human FGF-2/FGF basic/FGFb protein Met1-Ser155

Calculated MW 17.2 kDa
Observed MW 17 kDa
Accession BA G70135.1

**Bio-activity** Measured in a cell proliferation assay using BALB/c 3T3 cells. The ED<sub>50</sub> for this effect

is 0.75 ng/ml.

## **Properties**

**Purity** > 95 % as determined by reducing SDS-PAGE.

**Concentration** Subject to label value.

**Endotoxin**  $< 1.0 \text{ EU per } \mu\text{g}$  of the protein as determined by the LAL method.


Storage Storage Store at < -20°C, stable for 6 months. Please minimize freeze-thaw cycles.

**Shipping** This product is provided as liquid. It is shipped at frozen temperature with blue ice/gel


packs. Upon receipt, store it immediately at < - 20°C.

Formulation Supplied as a 0.2 μm filtered solution of 20mM Tris-HCl, 200mM NaCl, pH 7.5.

#### Data

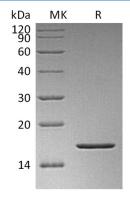

> 95 % as determined by reducing SDS-PAGE.

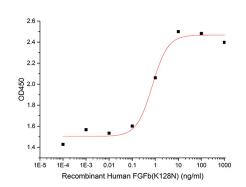

Measured in a cell proliferation assay using BALB/c 3T3 cells. The ED50 for this effect is 0.75 ng/ml.

## Background

Fibroblast growth factors (FGF) are a family of heparin-binding secreted proteins that stimulate cell proliferation and differentiation in a wide variety of tissues. FGFs play important roles in diverse biological functions both in vivo and in vitro; including mitogenesis; cellular migration; differentiation; angiogenesis; and wound healing. Human embryonic stem cell (hESC) cultures require FGF basic (also known as FGF-2 or bFGF) in cell culture media to remain in an undifferentiated and pluripotent state. Thermostable FGF basic was engineered for enhanced stability in culture media; without modification of its biological function. FGF basic is a required component of stem cell culture media for maintaining cells in an undifferentiated state. Because FGF basic is unstable; daily media changes are needed. The thermostable FGF basic that supports a 2-day media change schedule; so no media changes are required over a weeken d. This thermostable FGF basic was more stable than FGF basic in biochemical studies; and maintained cell growth; pluripotency and differentiation potential with a 2-day feeding schedule.

#### For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086 Fax: 1-832-243-6017