

A Reliable Research Partner in Life Science and Medicine

PE Anti-Mouse TER-119 Antibody[TER-119]

Catalog Number: E-AB-F1125D

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description

Reactivity Mouse Host Rat

Isotype Rat IgG2b, κ
Clone No. TER-119

Isotype Control PE Rat IgG2b, κ Isotype Control[LTF-2] [Product E-AB-F09842D]

Conjugation PE

Conjugation Information PE is designed to be excited by the Blue (488 nm), Green (532 nm) and Yellow-Green

(561 nm) lasers and detected using an optical filter centered near 575 nm (e.g., a 585/42

nm bandpass filter).

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer.


Applications Recommended usage

FCM Each lot of this antibody is quality control tested by flow cytometric analysis. The amount

of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for

individual use.

Data

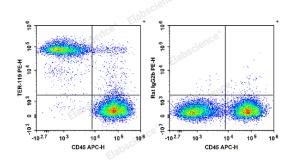

C57BL/6 murine bone marrow cells are stained with APC Anti-Mouse CD45 Antibody and PE Anti-Mouse TER-119 Antibody (Left). Bone marrow cells are stained with APC Anti-Mouse CD45 Antibody and PE Rat IgG2b, κ Isotype Control (Right).

Preparation & Storage

Storage Keep as concentrated solution.

This product can be stored at 2-8°C for 12 months. Please protected from prolonged

exposure to light and do not freeze.

Shipping Ice bag

Antigen Information

Alternate Names Ly-76;Lymphocyte antigen 76;TER119

Gene ID 104231

For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086 Fax: 1-832-243-6017

Web:www.elabscience.com Email:techsupport@elabscience.com

Elabscience Bionovation Inc.

A Reliable Research Partner in Life Science and Medicine

Background

Elabscience®

The TER-119 antigen is a 52 kD glycophorin A-associated protein, also known as Ly-76. TER-119 is an erythroid-specific antigen expressed on early proerythroblasts to mature erythrocytes, but not on erythroid colony-forming cells (BFU-E, blast-forming unit erythroid, or CFU-E, colony-forming unit erythroid).

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086
 Fax: 1-832-243-6017

 Web:www.elabscience.com
 Email:techsupport@elabscience.com