

Purified Anti-Human CD279/PD-1 Antibody[J116], Functional Grade

catalog number: E-AB-F11290

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description


Reactivity	Human
Immunogen	Recombinant Human CD279 protein
Host	Mouse
Isotype	Mouse IgG1, κ
Clone	J116
Purification	>98%, Protein A/G purified
Buffer	Sterile PBS, pH 7.2. < 1.0 EU per mg of the antibody as determined by the LAL method.

Applications


Recommended Dilution

FCM	$\leq 0.2 \mu\text{g}$ per million cells in 100 μL volume
Neut	Reported in the literature

Data

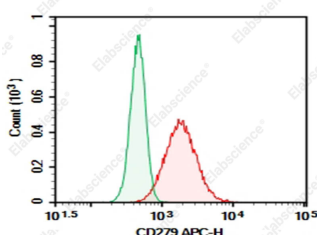

Molt-4 were activated for 4 hour with MIX, then stained with 0.2 μg Purified Anti-Human CD279 Antibody[J116], Functional Grade (Right) and 0.2 μg Mouse IgG1, κ Isotype Control(Left), followed by APC-conjugated Goat Anti-Mouse IgG Secondary Antibody.

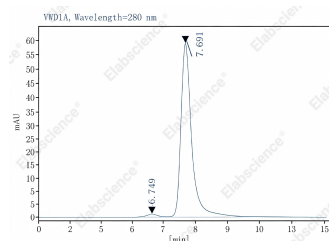

Molt-4 were activated for 4 hour with MIX, then stained with 0.2 μg Purified Anti-Human CD279 Antibody[J116], Functional Grade (Right) and 0.2 μg Mouse IgG1, κ Isotype Control(Left), followed by APC-conjugated Goat Anti-Mouse IgG Secondary Antibody.

Preparation & Storage

Storage	Store at 4°C valid for 12 months or -20°C valid for long term storage, avoid freeze / thaw cycles. This preparation contains no preservatives, thus it should be handled under aseptic conditions.
Shipping	Ice bag

Background

For Research Use Only

Toll-free: 1-888-852-8623
Web: www.elabscience.com

Tel: 1-832-243-6086
Email: techsupport@elabscience.com

Fax: 1-832-243-6017

Rev. V2.7

CD279, also known as programmed death-1 (PD-1), is a 50-55 kD glycoprotein belonging to the CD28 family of the Ig superfamily. PD-1 is expressed on activated splenic T and B cells and thymocytes. It is induced on activated myeloid cells as well. PD-1 is involved in lymphocyte clonal selection and peripheral tolerance through binding its ligands, B7-H1 (PD-L1) and B7-DC (PD-L2). It has been reported that PD-1 and PD-L1 interactions are critical to positive selection and play a role in shaping the T cell repertoire. PD-L1 negative costimulation is essential for prolonged survival of intratesticular islet allografts.

None (Azide-Free, Low Endotoxin) are perfectly suited to be used in culture or in vivo (for nonhuman studies) for functional assays blocking, neutralizing, activation or depletion where the presence of azide may damage cells or exogenous endotoxin may signal or activate cells.

Application References

Victor Tkachev, et al. J Immunol. 2015 Jun 15;194(12):5789-800.