
PerCP/Cyanine5.5 Anti-Human CD57 Antibody[HNK-1]

Catalog Number: E-AB-F1067J

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description	
Reactivity	Human
Host	Mouse
Isotype	Mouse IgM, к
Clone No.	HNK-1
Isotype Control	PerCP/Cyanine5.5 Mouse IgM, κ Isotype Control[MM-30] [Product E-AB-F09782J]
Conjugation	PerCP/Cyanine 5.5
Conjugation Information	PerCP/Cyanine5.5 is designed to be excited by the blue laser (488 nm) and detected using an optical filter centered near 675 nm (e.g., a 690/50 nm bandpass filter).
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Applications	Recommended usage
FCM	Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 µL of antibody per test (million cells in 100 µL staining volume or per 100 µL of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Data

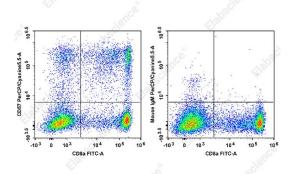

Staining of normal human peripheral blood cells with FITC Anti-Human CD8a Antibody[OKT-8] and PerCP/Cyanine5.5 Anti-Human CD57 Antibody[HNK-1] (left) or PerCP/Cyanine5.5 Mouse IgM, κ Isotype Control (right). Cells in the lymphocytes gate were used for analysis.

Keep as concentrated solution.
This product can be stored at 2-8°C for 12 months. Please protected from prolonged exposure to light and do not freeze.
lce bag
HNK-1;Leu-7;NK-1
Q9P2W7
e

For Research Use Only

Toll-free: 1-888-852-8623 Web:<u>w w w .elabscience.com</u>

Elabscience®

Gene ID Background

27087

CD57, also known as HNK-1, NK-1, and Leu-7 is a 100-115 kD oligosaccharide antigenic determinant expressed on a variety of proteins, lipids, and chondroitin sulfate proteoglycans. CD57 is expressed on a subset of peripheral blood lymphocytes, including NK cells and CD8+ T cells, and is also expressed on neural cells and striated muscle. CD57 is not expressed on red blood cells, granulocytes, monocytes, or platelets. While the function of CD57 is unknown, binding to L-selectin, P-selectin, and a fragment of laminin suggests that CD57 may be involved in cell-matrix interactions. CD57 is increased in some disease states associated with CD4/CD8 imbalances (AIDS, autoimmune disease, viral infections, and allograft transplants).

For Research Use Only