

A Reliable Research Partner in Life Science and Medicine

Elab Fluor® 647 Anti-Human CD45RO Antibody[UCHL1]

Catalog Number: E-AB-F1139M

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description

Reactivity Human Mouse Host

Isotype Mouse IgG2a, ĸ

Clone No. UCHL1

Isotype Control Elab Fluor[®] 647 Mouse IgG2a, κ Isotype Control[C1.18.4] [Product E-AB-F09802M]

Conjugation Elab Fluor®647

Conjugation Information Elab Fluor® 647 is designed to be excited by the Red laser (627-640 nm) and detected

using an optical filter centered near 670 nm (e.g., a 660/20 nm bandpass filter).


Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer. Storage Buffer

Applications Recommended usage

FCM

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 µL of antibody per test (million cells in 100 μL staining volume or per 100 μL of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Data

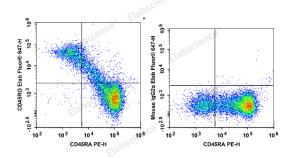

Human peripheral blood lymphocytes are stained with PE

Anti-Human CD45RA Antibody and Elab Fluor® 647 Anti-Human CD45RO Antibody (Left). Lymphocytes are stained

with PE Anti-Human CD45RA Antibody and Elab Fluor® 647 Mouse IgG2a, κ Isotype Control (Right).

Preparation & Storage

Storage Keep as concentrated solution.

This product can be stored at 2-8°C for 12 months. Please protected from prolonged

exposure to light and do not freeze.

Shipping Ice bag

Antigen Information

Gene ID 5788

For Research Use Only

Fax: 1-832-243-6017 Tel: 1-832-243-6086 Toll-free: 1-888-852-8623 Web:www.elabscience.com

Elabscience Bionovation Inc.



A Reliable Research Partner in Life Science and Medicine

Background

CD45RO is a 180 kD single chain membrane glycoprotein. It is a splice variant of tyrosine phosphatase CD45, lacking the A, B, and C determinants. The CD45RO isoform is expressed on activated and memory T cells, some B cell subsets, activated monocytes/macrophages, and granulocytes. CD45RO enhances both T cell receptor and B cell receptor signaling mediated activation. CD45 and its isoforms non-covalently associate with lymphocyte phosphatase-associated phosphoprotein (LPAP) on T and B lymphocytes. CD45 has been reported to be associated with several other cell surface antigens including CD1, CD2, CD3, and CD4. CD45 has also been reported to bind galectin-1 and CD22. CD45 isoform expression can change in response to cytokines.

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086
 Fax: 1-832-243-6017

 Web:www.elabscience.com
 Email:techsupport@elabscience.com