

icella



# Ham's F-12 Complete Medium (20)

Cat. No.: PM150810C

Size: 125mL×4

#### **General Information**

**Product Form** Liquid Concentration

Ready-to-use Components Ham's F-12[PM150810]+Nutrients+2 mM L-Alanyl-L-Glutamine[PB180419]

**Bacterial detection** Negative Fungal detection Negative Mycoplasmal detection Negative **Endotoxin level** < 3 EU/mL Shipping Ice bag

2-8°C, Shading Light Storage

**Expiry date** 3 months

### **Product Introduction**

Ham's F-12 Nutrient Mixture was designed by Ham in 1969 based on Ham's F-10 nutrient mixture, which was originally used for serum-free culture of CHO cells. Ham's F-12 is often used as the basic culture medium for serum free culture. At low serum levels, it is especially suitable for single cell culture and clonal culture. After adding serum, it is also widely used in the culture of cancer cells and primary cells, such as rat hepatocytes, rat prostate epithelial cells, chondrocytes, rat myoblasts, and chicken embryonic cells. In addition, the DMEM/F12 medium obtained by mixing Ham's F-12 in equal volume with DMEM is more nutritious and widely used.

## **Guidelines for use**

- 1. Pricella's cell culture media undergoes strictquality control to ensure sterility, but may get contaminated during use. Follow these guidelines for sterile handling toavoid contamination.
- 2. Always wipe your gloved hands and work area with 70% ethanol.
- 3. Wipe the outside of the containers, flasks, plates, and dishes with 70% ethanol before placing them in the cell culture hood.
- 4. Use sterile pipette tips and pipettes to work with liquids, and use each pipette tip only once to avoid cross contamination.Do not unwrap sterile pipettes until they are ready to beused. Keep pipettes and tips within the clean work area.
- 5. Do not talk while performing sterile procedures and perform your cell culture as efficiently and carefully as possible tominimize contamination.

## **Quality control**

Standard evaluations for cell culture media are pH, osmolality, endotoxins and sterility testing for liquid products, by Elabscience cell growth experiments.


**Notes** 

Tel:1-832-243-6086 Fax:1-832-243-6017 Toll-free:1-888-852-8623 Email: techsupport@pri-cella.com Web: www.pri-cella.com





- 1. This product is for research use only.
- 2. This product is sterilized by 0.1  $\mu m$  filtration.
- 3. It is necessary to pay attention to the aseptic operation and avoid the contamination during the culture.
- 4. It is not suitable for long time storage at room temperature.
- 5. This product is a ready-to-use medium. If there is no special need, don't add serum, penicillin and streptomycin. It can be used directly.

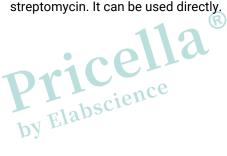







Tel:1-832-243-6086 Fax:1-832-243-6017 Toll-free:1-888-852-8623 Web: www.pri-cella.com