

A Reliable Research Partner in Life Science and Medicine

Elab Fluor® 700 Anti-Human CD329 Antibody[K8]

Catalog Number: AN00319M1

Note: Centrifuge before opening to ensure complete recovery of vial contents.

_	e				40		
		•	^		TI.	$\boldsymbol{\cap}$	
	4-	-	•	w	ш		

Reactivity Human Host Mouse

Isotype Mouse IgG1, κ

Clone No. K8

Isotype Control Elab Fluor® 700 Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792M1]

Conjugation Elab Fluor® 700

Conjugation Information Elab Fluor® 700 is designed to be excited by the Red laser (627-640 nm) and detected

using an optical filter centered near 719 nm (e.g., a 725/40 nm bandpass filter).


Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% sodium azide and 1% BSA.

Applications Recommended usage

FCM

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Data

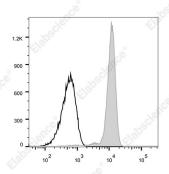

Staining of normal human peripheral blood cells with Elab

Fluor[®] 700 Anti-Human CD329 Antibody[K8](filled gray histogram) or Elab Fluor[®] 700 Mouse IgG1, κ Isotype Control(empty black histogram). Cells in the granulocytes gate were used for analysis.

Preparation & Storage

Storage Keep as concentrated solution.

This product can be stored at 2-8°C for 12 months. Please protected from prolonged

exposure to light and do not freeze.

Web: www.elabscience.cn

Shipping Ice bag

Antigen Information

Alternate Names Siglec-9; Sialic acid-binding Ig-like lectin-9

 Uniprot ID
 Q9Y336

 Gene ID
 27180

For Research Use Only



Elabscience Biotechnology Co., Ltd.

A Reliable Research Partner in Life Science and Medicine

Background

Siglecs are cell surface receptors belonging to the immunoglobulin superfamily that recognize sugar antigens. The extracellular domain of siglec-9 contains an IgV region, which binds sialic acid, followed by two IgC regions. Siglec 9 and siglec 6-8,10-12 are CD33 (siglec 3) like siglecs, which have two ITIMs in the cytoplasmic tails, suggesting their functional involvement in signal transduction. It is highly expressed on neutrophils and monocytes, and at lower levels on the subpopulations of T and B lymphocytes and NK cells. Siglec-9 plays a role in negative regulation of T cell activation, and it also affects neutrophil apoptosis.

Rev. V1.4

Web: www.elabscience.cn