

(FOR RESEARCH USE ONLY. DO NOT USE IT IN CLINICAL DIAGNOSIS !)

**Catalog No: E-BC-K176-M**

**Specification: 48T(32 samples)/96T(80 samples)**

**Measuring instrument: Microplate reader(580-590 nm)**

**Detection range: 0.70-80  $\mu\text{mol/L}$**

## **Elabsience® Lipid Peroxide (LPO) Colorimetric Assay Kit**

This manual must be read attentively and completely before using this product.

If you have any problem, please contact our Technical Service Center for help:

Toll-free: 1-888-852-8623

Tell: 1-832-243-6086

Fax: 1-832-243-6017

Email: [techsupport@elabsience.com](mailto:techsupport@elabsience.com)

Website: [www.elabsience.com](http://www.elabsience.com)

Please kindly provide us the lot number (on the outside of the box) of the kit for more efficient service.

## Table of contents

<b>Assay summary .....</b>	<b>3</b>
<b>Intended use .....</b>	<b>4</b>
<b>Detection principle .....</b>	<b>4</b>
<b>Kit components &amp; storage .....</b>	<b>4</b>
<b>Materials prepared by users .....</b>	<b>5</b>
<b>Reagent preparation .....</b>	<b>5</b>
<b>Sample preparation .....</b>	<b>6</b>
<b>The key points of the assay .....</b>	<b>7</b>
<b>Operating steps .....</b>	<b>7</b>
<b>Calculation .....</b>	<b>8</b>
<b>Appendix I Performance Characteristics .....</b>	<b>9</b>
<b>Appendix II Example Analysis .....</b>	<b>11</b>
<b>Appendix III Publications .....</b>	<b>12</b>
<b>Statement .....</b>	<b>13</b>

**Assay summary**



## Intended use

This kit can be used to measure LPO content in serum, plasma, urine and tissue amples.

## Detection principle

With 45°C incubation for 60 min, one molecule of LPO react with two molecule of chromogenic reagent, to produce a stable chromophore which have the maximum absorption peak at 586 nm. The content of LPO in samples can be calculated by standard curve or calculation formula.

## Kit components & storage

Item	Component	Size 1(48 T)	Size 2(96 T)	Storage
Reagent 1	Substrate Stock Solution	30 mL × 1 vial	60 mL × 1 vial	2-8°C, 12 months, shading light
Reagent 2	Diluent	20 mL × 1 vial	20 mL × 1 vial	2-8°C, 12 months, shading light
Reagent 3	Acid Reagent	10 mL × 1 vial	20 mL × 1 vial	2-8°C, 12 months, shading light
Reagent 4	100 µmol/L Standard	6 mL ×1 vial	6 mL ×1 vial	2-8°C, 12 months, shading light
	Microplate	48 wells	96 wells	No requirement
	Plate Sealer	2 pieces		
	Sample Layout Sheet	1 piece		

Note: The reagents must be stored strictly according to the preservation conditions in the above table. The reagents in different kits cannot be mixed with each other. For a small volume of reagents, please centrifuge before use, so as not to obtain sufficient amount of reagents.

## Materials prepared by users

### Instruments:

Microplate reader (580–590 nm, optimum wavelength: 586 nm),

Micropipettor, Centrifuge, Incubator, Vortex mixer.

### Reagents:

Anhydrous ethanol, PBS (0.01 M, pH 7.4)

## Reagent preparation

① Equilibrate all the reagents to room temperature before use.

② The preparation of chromogenic working solution:

For each well, prepare 650  $\mu\text{L}$  of chromogenic working solution (mix well 487.5  $\mu\text{L}$  of substrate stock solution and 162.5  $\mu\text{L}$  of diluent). The chromogenic working solution should be prepared on spot.

③ The preparation of standard curve:

Always prepare a fresh set of standards. Discard working standard dilutions after use.

Dilute 100  $\mu\text{mol/L}$  standard solution with absolute ethanol to a serial concentration. The recommended dilution gradient is as follows: 0, 5, 10, 20, 30, 40, 50, 80  $\mu\text{mol/L}$ . Reference is as follows:

Item	①	②	③	④	⑤	⑥	⑦	⑧
<b>Concentration (<math>\mu\text{mol/L}</math>)</b>	<b>0</b>	<b>5</b>	<b>10</b>	<b>20</b>	<b>30</b>	<b>40</b>	<b>50</b>	<b>80</b>
<b>100 <math>\mu\text{mol/L}</math> standard (<math>\mu\text{L}</math>)</b>	0	50	100	200	300	400	500	800
<b>Absolute ethanol (<math>\mu\text{L}</math>)</b>	1000	950	900	800	700	600	500	200

## Sample preparation

### ① Sample preparation:

**Serum and plasma:** detect directly. If not detected on the same day, the serum or plasma can be stored at  $-80^{\circ}\text{C}$  for a month.

**Urine:** Collect fresh urine and centrifuge at  $10000\times g$  for 15 min at  $4^{\circ}\text{C}$ . Take the supernatant to preserve it on ice for detection. If not detected on the same day, the urine can be stored at  $-80^{\circ}\text{C}$  for a month.

### Tissue samples:

- ① Harvest the amount of tissue needed for each assay (initial recommendation 40 mg).
- ② Wash tissue in cold PBS (0.01 M, pH 7.4).
- ③ Homogenize 40 mg tissue in 360  $\mu\text{L}$  PBS (0.01 M, pH 7.4) with a dounce homogenizer at  $4^{\circ}\text{C}$ .
- ④ Centrifuge at  $10000\times g$  for 10 min to remove insoluble material. Collect supernatant and keep it on ice for detection.
- ⑤ Meanwhile, determine the protein concentration of supernatant (E-BC-K318-M).

### ② Dilution of sample

The recommended dilution factor for different samples is as follows (for reference only):

Sample type	Dilution factor
Human serum	1
Human plasma	1
Urine	1
10% Rat kidney tissue homogenization	1
10% Rat liver tissue homogenization	1

Note: The diluent is PBS (0.01 M, pH 7.4). For the dilution of other sample types, please do pretest to confirm the dilution factor.

## The key points of the assay

- ① The EP tube needs to be sealed to avoid leakage.
- ② The supernatant added to 96-well microplate must be clarified, otherwise centrifuge again.
- ③ Please carry out the experiment in the fume hood and wear disposable gloves, prevent the reagents splashing dashed into eyes and skin.

## Operating steps

- ① Standard well: add 200  $\mu$ L of standard solution with different concentration into the 1.5 mL EP tube.  
Sample well: add 200  $\mu$ L of sample into the 1.5 mL EP tube.
- ② Add 650  $\mu$ L of chromogenic working solution, cover the caps and mix fully.
- ③ Add 150  $\mu$ L of acid reagent, cover the caps and mix fully.
- ④ Incubate at 45°C for 60 min. Cool to room temperature with running water.
- ⑤ Centrifuge at 1100 $\times$ g for 10 min. Take 200  $\mu$ L of supernatant to the microplate, measure the OD values of each well at 586 nm with microplate reader. (Avoid the bubbles generated when adding the liquid to microplate, otherwise OD value will be affected.)

## Calculation

### The standard curve:

1. Average the duplicate reading for each standard.
2. Subtract the mean OD value of the blank (Standard # ① ) from all standard readings. This is the absolved OD value.
3. Plot the standard curve by using absolved OD value of standard and correspondent concentration as y-axis and x-axis respectively. Create the standard curve ( $y = ax + b$ ) with graph software (or EXCEL).

### The sample:

#### 1. Serum (plasma) and other liquid sample:

$$\text{LPO } (\mu\text{mol/L}) = (\Delta A_{586} - b) \div a \times f$$

#### 2. Tissue sample (Calculated by tissue protein):

$$\text{LPO } (\mu\text{mol/gprot}) = (\Delta A_{586} - b) \div a \times f \div C_{pr}$$

#### 3. Tissue sample (Calculated by tissue wet weight):

$$\text{LPO } (\mu\text{mol/kg wet weight}) = (\Delta A_{586} - b) \div a \times f \div m \times V$$

### [Note]

$\Delta A_{586}$ : Absolute OD ( $OD_{\text{Sample}} - OD_{\text{Blank}}$ ).

f: Dilution factor of sample before test.

$C_{pr}$ : Concentration of protein in sample (gprot/L).

m: The wet weight of tissue, g.

V: The volume of homogenate of tissue sample, mL.



## Appendix I Performance Characteristics

### 1. Parameter:

#### Intra-assay Precision

Three human serum samples were assayed in replicates of 20 to determine precision within an assay. (CV = Coefficient of Variation)

Parameters	Sample 1	Sample 2	Sample 3
Mean ( $\mu\text{mol/L}$ )	5.60	34.50	66.00
%CV	3.5	2.9	2.9

#### Inter-assay Precision

Three human serum samples were assayed 20 times in duplicate by three operators to determine precision between assays.

Parameters	Sample 1	Sample 2	Sample 3
Mean ( $\mu\text{mol/L}$ )	5.60	34.50	66.00
%CV	3.2	3.5	3.8

#### Recovery

Take three samples of high concentration, middle concentration and low concentration to test the samples of each concentration for 6 times parallelly to get the average recovery rate of 99%.

	Standard 1	Standard 2	Standard 3
Expected Conc. ( $\mu\text{mol/L}$ )	7.5	26.4	48.5
Observed Conc. ( $\mu\text{mol/L}$ )	7.6	26.1	47.0
Recovery rate (%)	101	99	97


#### Sensitivity

The analytical sensitivity of the assay is  $0.70 \mu\text{mol/L}$ . This was determined by adding two standard deviations to the mean O.D. obtained when the zero standard was assayed 20 times, and calculating the corresponding concentration.

## 2. Standard curve

As the OD value of the standard curve may vary according to the conditions of the actual assay performance (e.g. operator, pipetting technique or temperature effects), so the standard curve and data are provided as below for reference only.

Concentration ( $\mu\text{mol/L}$ )	0	5	10	20	30	40	50	80
Average OD	0.079	0.140	0.236	0.389	0.563	0.722	0.903	1.357
Absoluted OD	0	0.061	0.157	0.31	0.484	0.643	0.824	1.278

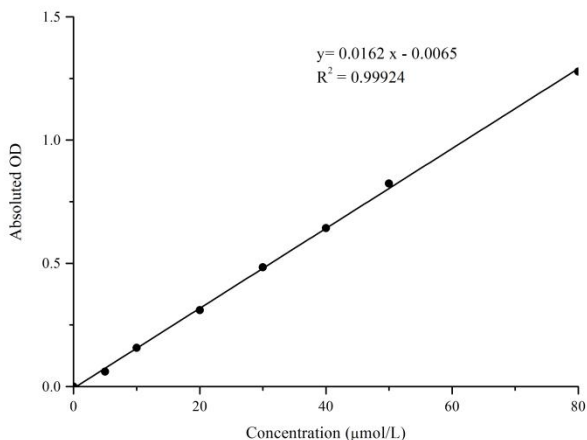

## Appendix Π Example Analysis


### Example analysis :

Take 200  $\mu\text{L}$  of human serum, carry the assay according to the operation steps. The results are as follows:

standard curve:  $y = 0.0162x - 0.0065$ , the average OD value of the sample is 0.219, the average OD value of the blank is 0.073, and the calculation result is:

$$\text{LPO } (\mu\text{mol/L}) = (0.219 - 0.073 + 0.0065) \div 0.0162 \times 1 = 9.41 \text{ } (\mu\text{mol/L})$$

Detect human serum, human urine, 10% rat liver tissue homogenate (the concentration of protein is 13.10 gprot/L), 10% rat kidney tissue homogenate (the concentration of protein is 9.26 gprot/L) according to the protocol, the result is as follows:

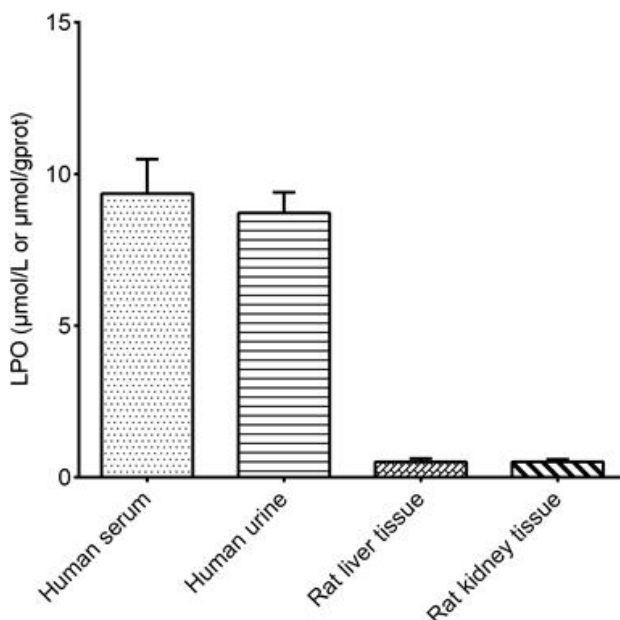

## Appendix III Publications

1. Zhang L , Shi W Y ,Jia-Ying XuYan LiuShi-Jia WangJia-Yang ZhengYun-Hong LiLin-Xi YuanLi-Qiang Qin.Protective effects and mechanism of chemical- and plant-based selenocystine against cadmium-induced liver damage[J].Journal of Hazardous Materials, 2024, 468(Apr.15):133812.1-133812.12.DOI:10.1016/j.jhazmat.2024.133812.
2. Wu Z , Chen K , Mo W ,et al.Multimodal enhancement of ferroptosis for synergistic cascade colorectal cancer therapy[J].Chemical Engineering Journal, 2024, 498(000):15.DOI:10.1016/j.cej.2024.155048.
3. Wang D , Zhang H , Liao X ,et al.Oral administration of Robinia pseudoacacia L. flower exosome-like nanoparticles attenuates gastric and small intestinal mucosal ferroptosis caused by hypoxia through inhibiting HIF-1 $\alpha$ - and HIF-2 $\alpha$ -mediated lipid peroxidation[J].Journal of Nanobiotechnology, 2024, 22(1):1-19.DOI:10.1186/s12951-024-02663-6.
4. Leng B , Deng L , Tan J ,et al.Targeting the Na<sup>+</sup>/K<sup>+</sup> ATPase DR-region with DR-Ab improves doxorubicin-induced cardiotoxicity[J].Free Radical Biology & Medicine, 2023, 204(000):16.DOI:10.1016/j.freeradbiomed.2023.04.008.
5. Xu J , Ding Y , Shi C ,et al.Identification of Cathepsin B as a Therapeutic Target for Ferroptosis of Macrophage after Spinal Cord Injury[J].Aging & Disease, 2024, 15(1).DOI:10.14336/AD.2023.0509.
6. Yang L , Zhang Z , Yao X ,et al.HNRNPL facilitates ferroptosis in hepatocellular carcinoma cells by promoting S100A9 expression[J].Translational Oncology, 2024, 43.DOI:10.1016/j.tranon.2024.101908.

## Statement

1. This assay kit is for Research Use Only. We will not response for any arising problems or legal responsibilities causing by using the kit for clinical diagnosis or other purpose.
2. Please read the instructions carefully and adjust the instruments before the experiments. Please follow the instructions strictly during the experiments.
3. Protection methods must be taken by wearing lab coat and latex gloves.
4. If the concentration of substance is not within the detection range exactly, an extra dilution or concentration should be taken for the sample.
5. It is recommended to take a pre-test if your sample is not listed in the instruction book.
6. The experimental results are closely related to the situation of reagents, operations, environment and so on. Elabscience will guarantee the quality of the kits only, and NOT be responsible for the sample consumption caused by using the assay kits. It is better to calculate the possible usage of sample and reserve sufficient samples before use.





