

A Reliable Research Partner in Life Science and Medicine

PerCP Anti-Human CD45RA Antibody[HI100]

Catalog Number: E-AB-F1052F

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description

Reactivity Human Host Mouse

Isotype Mouse IgG2b, κ

Clone No. HI100

Isotype Control [Product E-AB-F09812F]

Conjugation PerCP

Conjugation Information PerCP is designed to be excited by the blue laser (488 nm) and detected using an optical

filter centered near 675 nm (e.g., a 690/50 nm bandpass filter).


Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% sodium azide and 1% BSA.

Applications Recommended usage

FCM

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Data

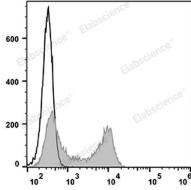

Human peripheral blood lymphocytes are stained with PerCP Anti-Human CD45RA Antibody (filled gray histogram). Unstained lymphocytes (empty black histogram) are used as control.

Preparation & Storage

Storage Keep as concentrated solution.

This product can be stored at 2-8°C for 12 months. Please protected from prolonged

exposure to light and do not freeze.

Shipping Ice bag

Antigen Information

Alternate Names CD45;L-CA;Ly-5;Ptprc;Receptor-type tyrosine-protein phosphatase C;T200

Web: www.elabscience.cn

 Uniprot ID
 P08575

 Gene ID
 5788

For Research Use Only

Elabscience®

Elabscience Biotechnology Co., Ltd.

A Reliable Research Partner in Life Science and Medicine

Background

CD45RA is a 205-220 kD single chain type I glycoprotein. It is an exon 4 splice variant of the tyrosine phosphatase CD45. The CD45RA isoform is expressed on resting/naïve T cells, medullary thymocytes, B cells and monocytes. CD45RA enhances both T cell receptor and B cell receptor signaling. CD45 non-covalently associates with lymphocyte phosphatase-associated phosphoprotein (LPAP) on T and B lymphocytes. CD45 has been reported to be associated with several other cell surface antigens including CD1, CD2, CD3, and CD4. CD45 has also been reported to bind galectin-1. CD45 isoform expression can change in response to cytokines.

Web: www.elabscience.cn