



McCoy's 5A (with HEPES) Complete Medium (20)

Cat. No.: PM150711C

Size: 125mL×4

General Information

Product Form Liquid Concentration Ready-to-use

Components McCoy's 5A, with HEPES[PM150711]+20% Nutrient+1% Supplement

Bacterial detection Negative Fungal detection Negative Mycoplasmal detection Negative **Endotoxin level** < 3 EU/mL Shipping Ice bag

2-8°C, Shading Light Storage

Expiry date 6 months

Product Introduction

icella McCoy's 5A medium was initially used specifically for the culture of Novikoff Hepatoma cells. Unlike other medium, McCoy's 5A medium contains reduced glutathione, bacterial peptone and high concentrations of glucose. McCoy's 5A is widely used in the culture of primary cells, such as bone marrow, skin, gingiva, kidney, spleen, lung, rat embryo, omentum and so on. In addition, McCoy's 5A medium is also used for tissue biopsy culture, the establishment of cell line, and the culture of some lymphocytes and more difficult cells.

Guidelines for use

- 1. Pricella's cell culture media undergoes strictguality control to ensure sterility, but may get contaminated during use. Follow these guidelines for sterile handling toavoid contamination.
- 2. Always wipe your gloved hands and work area with 70% ethanol.
- 3. Wipe the outside of the containers, flasks, plates, and dishes with 70% ethanol before placing them in the cell culture hood.
- 4. Use sterile pipette tips and pipettes to work with liquids, and use each pipette tip only once to avoid cross contamination. Do not unwrap sterile pipettes until they are ready to beused. Keep pipettes and tips within the clean work area.
- 5. Do not talk while performing sterile procedures and perform your cell culture as efficiently and carefully as possible tominimize contamination.

Quality control

PT1CEIIIty te Standard evaluations for cell culture media are pH, osmolality, endotoxins and sterility testing for liquid products, cell growth experiments.


Notes

Tel:1-832-243-6086 Fax:1-832-243-6017 Toll-free:1-888-852-8623 Email: techsupport@pri-cella.com Web: www.pri-cella.com





- 1. This product is for research use only.
- 2. This product is sterilized by 0.1 μm filtration.
- 3. It is necessary to pay attention to the aseptic operation and avoid the contamination during the culture.
- 4. It is not suitable for long time storage at room temperature.
- 5. This product is a ready-to-use medium. If there is no special need, don't add serum, penicillin and streptomycin. It can be used directly.

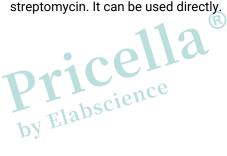







Tel:1-832-243-6086 Fax:1-832-243-6017 Toll-free:1-888-852-8623 Web: www.pri-cella.com