A Reliable Research Partner in Life Science and Medicine

# Recombinant Mouse E-Selectin/SELE Protein (Fc Tag)

Catalog Number: PKSM040510

Note: Centrifuge before opening to ensure complete recovery of vial contents.

#### **Description**

**Species** Mouse

Source HEK293 Cells-derived Mouse E-Selectin/SELE protein Met1-Pro564, with an C-

terminal hFc

Calculated MW 85.5 kDa Observed MW 101 kDa Accession NP 035475.1

Measured by the ability of the immobilized protein to support the adhesion of U937 **Bio-activity** 

> human histiocytic lymphoma cells. When cells are added to coated plates (10 µg/ml, 100 μl/well), approximately 50% will adhere specifically after 60 minutes at 37 °C.

### **Properties**

Purity > 90 % as determined by reducing SDS-PAGE.

Endotoxin < 1.0 EU per µg of the protein as determined by the LAL method.

Generally, lyophilized proteins are stable for up to 12 months when stored at -20 to -80 Storage

°C. Reconstituted protein solution can be stored at 4-8°C for 2-7 days. Aliquots of

reconstituted samples are stable at < -20°C for 3 months.

This product is provided as lyophilized powder which is shipped with ice packs. Shipping

Lyophilized from sterile PBS, pH 7.4 Formulation


Normally 5% - 8% trehalose, mannitol and 0.01% Tween 80 are added as protectants

before lyophilization.

Please refer to the specific buffer information in the printed manual.

Reconstitution Please refer to the printed manual for detailed information.

#### Data

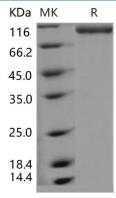

> 90 % as determined by reducing SDS-PAGE.

## Background

#### Elabscience Bionovation Inc.



A Reliable Research Partner in Life Science and Medicine

E-selectin, also known as endothelial leukocyte adhesion molecule-1 (ELAM-1) and CD62E, is an inducible adhesion molecule that is expressed on the surfaces of stimulated vascular endothelial cells and is sometimes involved in cancer cell metastasis. As a member of the LEC-CAM or selectin family, E-selectin recognises and binds to sialylated carbohydrates including members of the Lewis X and Lewis A families found on monocytes, granulocytes, and T-lymphocytes. E-selectin supports rolling and stable arrest of leukocytes on activated vascular endothelium, and furthermore, it was indicated that it can also transduce an activating stimulus via the MAPK cascade into the endothelial cell during leukocyte adhesion. E-selectin regulates adhesive interactions between certain blood cells and endothelium. The soluble form of E selectin (sE-selectin) is a marker of endothelial activation, and has a potential role in the pathogenesis of cardiovascular disease as raised levels have been found in hypertension, diabetes and hyperlipidemia, although its association in established atherosclerosis disease and its value as a prognostic factor is more controversial. soluble E-selectin is inversely associated with the muscular component of the left ventricle, thereby suggesting that the lack of such a reparative factor may be associated with cardiac remodeling in end-stage renal disease (ESRD) patients. In addition, this adhesion molecule appears to be involved in the pathogenesis of atherosclerosis.

Fax: 1-832-243-6017